Mode d'emploi

Instruction manual

VEINOPUS

Model 1.2

SOMMAIRE

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www.adremtechnology.com

SUMMARY

• SYMBOLS	12
• FUNCTION AND INTENDED APPLICATION	12
• INDICATIONS AND BENEFITS	12
• CONTRA-INDICATIONS AND WARNINGS	13
• SAFETY	13
• HOW TO USE VEINOPLUS	13
• TECHNICAL DESCRIPTION	15
• ACCESSORIES	16
• RISKS ASSOCIATED WITH DISPOSAL OF ACCESSORIES OR PRODUCT	17
• ENVIRONMENTAL CONDITIONS FOR USE, STORAGE AND TRANSPORT	17
• SERVICE	18
• WARRANTY	18
• PICTURE	19

SYMBOLS

Read the instructions before using the device

A device type BF (IEC standard N° 601-1) with an internal power supply and intended for continuous use with IPX0 degree of protection against water penetration

Cable outlet indication (4 pole RJ 45)

Adjustment of intensity

ON/OFF switch

FUNCTION AND INTENDED APPLICATION

The **function** of VEINOPLUS is to electrically stimulate muscles via motor nerves. Its **intended application** is to cause contractions of calf muscles in legs. These calf contractions, which feel like pulsations, compress deep veins and result in pumping of venous blood against gravity.

Who can benefit from VEINOPLUS ?

People suffering from swelling and pain sensation in legs due to:

- Prolonged sitting or standing (office, air travel)
- Lack of physical activity
- Pregnancy
- Hormone therapy
- Body overheating
- Obesity

INDICATIONS AND BENEFITS

VEINOPLUS improves venous blood return, which results in:

- Control of pain sensations in legs
- Reduction of swelling (œdemas) of feet and ankles
- Prevention of Deep Venous Thrombosis (DVT)
- Prevention of venous stasis

CONTRA-INDICATIONS AND WARNINGS

Contra-indications:

- People with **implanted cardiac pacemakers**.
- Patients undergoing anti-cancer therapy.

Warnings:

If you are suffering from one of the following conditions or if in doubt, consult a physician before using VEINOPLUS:

- Acute venous thrombosis
- Tendency to bleeding induced by anti - clotting treatment
- Sciatica or any other neuralgia
- Inflammation, infection, skin lesions, or severe varicose veins of the calves: in these instances you may apply the electrodes next to, but not onto the affected area

Precautions and limitations:

- The use of VEINOPLUS is limited to the legs as well as arms in certain cases
- It is not intended for application on head, chest, abdomen or neck
- Use of VEINOPLUS near a heart region may increase the risk of ventricular fibrillation
- Keep out of reach of children

SAFETY

The VEINOPLUS produces electrical impulses and delivers energy to the surface of the skin by means of its electrodes, which cause local muscular contractions. The impulses from VEINOPLUS carry very low quantities of electrical energy. Their energies are below any international standard safety limits. The use of unauthorized accessories or materials may render VEINOPLUS ineffective, but it will not degrade its safety.

HOW TO USE VEINOPLUS (see picture on page 19)


- 1 **Connect the cable** to the outlet located on the top of the unit (7)
Connect the electrodes to the electrode end (the one with pins) of the cable
- 2 **Stick each of the electrodes** to the central part of the calf muscle on the back of legs

- 3 **Make yourself comfortable:** the treatment lasts 20 minutes. The ideal positions are: seated with the legs slightly extended or lying down
- 4 **Switch ON the stimulator** by pressing the ON/OFF button (1)
- 5 **The intensity may be adjusted** using the two push buttons (2). To increase intensity, press the button marked **+**, to decrease it, press the button marked **-**

You should increase the intensity until you feel light but deep muscular pulsations in your calves. They may vary between one calf and the other. The intensity should be adjusted in relation to the more sensitive of the two calves. You can decrease the intensity, if the pulsations become too strong.

- 6 **Duration:** 20 minutes, during which the rate of pulsations changes every 5 minutes. The VEINOPLUS switches OFF automatically after 20 minutes
- 7 **When to use :** daily or as needed (no limits, no risk of overuse). Especially recommended during and up to 2 weeks after long-haul flights
- 8 **OFF:** If for any reason you wish to switch OFF VEINOPLUS, push once the ON/OFF button (1)
- 9 **The effectiveness of VEINOPLUS depends on:**
 - The correct positioning and using electrodes marked VEINOPLUS PACK.
 - The correct adjustment of the intensity of stimulation, up to reaching light but deep pulsations of calf muscles.
 - If your skin is too dry, wet it with tap water, to further improve the effectiveness of stimulation.

The VEINOPLUS treatment may be even more effective when wearing knee high stockings or support socks during electrical stimulation. The electrodes must be placed directly on to the skin, and not on the stockings.

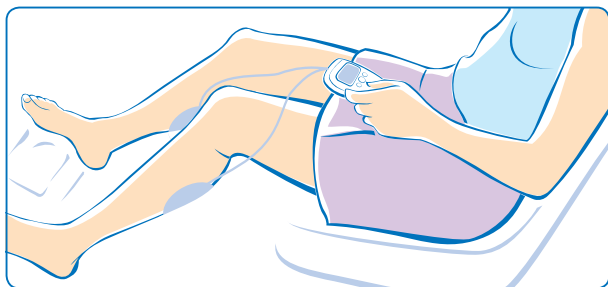TECHNICAL DESCRIPTION

VEINOPLUS is a battery powered electrical muscle stimulator that delivers specific low frequency and low voltage currents. These currents are a highly developed synthesis of impulses. They provide safe and effective muscle pulsations in legs to facilitate venous blood return against gravity. VEINOPLUS is powered by a single 9V alkaline battery. Input power does not exceed 0.3W. **All parameters of the output impulses do not change throughout the life of the battery.**

Features (see picture on page 19):

- One push button for ON / OFF (1)
- Two push buttons for intensity control (2)
- An LCD screen displaying the intensity of stimulation (3), - A vertical bar flashes each time the muscle pumping occurs-
- Remaining time left to automatic stop (4),
- A battery icon on the screen flashes when the battery needs to be replaced (5)
- Counter of utilisations (6)
- One single insulated stimulation channel (7)
- Two electrodes (8)

An internal timer switches the device OFF after 20 minutes

Output parameters of VEINOPLUS:

- The low frequency current is modulated in a special pattern from 1 Hz to 250 Hz
- The waveform of the current impulses during use is quasi-symmetrical biphasic AC.
- The voltage waveform of the impulse is rectangular and monopolar.
- The time duration of output impulses ranges from 25 to 240 microseconds.
- During impulses, the instantaneous output voltage reaches a peak value of $45V_p \pm 10\%$, when loaded with impedance in a range from 250 ohms to 10'000 ohms

Maximum ratings of the VEINOPLUS output when loaded with 500 Ω resistance:

- The effective value of AC output voltage is less than 2V r.m.s over 5 seconds.
- The effective value of AC current intensity does not exceed 4 mA r.m.s. over 5 seconds and it may reach a peak value of $90\text{ mA} \pm 10\%$ during impulses.

- The maximum r.m.s. value of current density on recommended electrodes does not exceed 0.2 mA/cm²
- The maximum AC output power of VEINOPLUS is limited to 0.01 W r.m.s.
- The maximum DC output power of VEINOPLUS is limited to 0.0001 W.

Rates of muscle pulsations during use:

The rates of stimulated muscle pulsations are similar to heart beats. They are programmed to accelerate during each use from 60 BPM (Beats Per Minute) for the first 5 minutes of treatment, to 75 BPM for the following 5 minutes, then to 90 BPM for the next 5 minutes and finally to 105 BPM for the last 5 minutes of use session. These pulsation rates correspond to frequencies in a range from 1 Hz to 1.75 Hz.

ACCESSORIES

VEINOPLUS is supplied with all the accessories required for an easy, safe and effective treatment. These include **two VEINOPLUS PACK self-adhesive skin electrodes, one cable, one 9V alkaline battery and one instruction manual.** In compliance with transport regulations in some countries, the battery may not be delivered with the device.

Electrodes and how to use them:

The VEINOPLUS PACK electrodes are reusable. It is recommended that the same set of electrodes is used by the same person. Used electrodes should be replaced as soon as they no longer adhere or become dirty. Specific instructions for use and maintenance of the electrodes are included in their package.

Recommended electrodes:

VEINO PLUS has been clinically tested, and its effectiveness is optimal when VEINOPLUS PACK 8X13 electrodes are used. For effective and safe treatments it is recommended to use only the original electrodes branded VEINOPLUS PACK. No other brands or types of electrodes are recommended for use with VEINOPLUS. The replacement electrodes when purchased in pharmacies or elsewhere should be in original sealed packages VEINOPLUS PACK. If, by mistake any other types of electrodes are used with the VEINOPLUS, its effectiveness can be severely reduced but its safety will not be affected. This is due to inherent design of electronic circuitry inside the VEINOPLUS, which automatically limits the effective current density on electrode contact if smaller or larger size auto-adhesive electrodes are used by mistake. Metal plates for skin contacts instead of pre-gelled VEINOPLUS PACK electrodes can be used only by, and at the sole responsibility of health professionals.

The battery and its replacement:

When the 9v battery reaches its end-of-life, (voltage < 6V), the battery icon on the screen begins to flash. You should change the battery either at the end of the session or when the VEINOPLUS switches OFF.

To replace the battery no tools are needed, just open the cover on the back of the VEINOPLUS and insert a new battery observing its polarity. It is recommended to use 9V alkaline batteries of any type listed below: 6LR61, 6LF22, AM6, MN1604 or GP 1604. If other batteries than those recommended are used, the VEINOPLUS may cease functioning effectively or prematurely, but its safety will not be affected.

The electrode cable and its replacement:

The cable supplied with VEINOPLUS can be replaced, if needed, through your local VEINOPLUS distributor.

RISKS ASSOCIATED WITH DISPOSAL OF ACCESSORIES OR PRODUCT

There are no risks associated with handling and disposal of waste created by VEINOPLUS product or its accessories, except for the batteries. The following precautions relate to batteries:

- The batteries must not be exposed to temperatures above 65° C
- When your device is not used for a long period of time, remove the battery
- The alkaline battery used in VEINOPLUS must not be recharged as it may explode
- The battery must not be disposed into fire as it may explode
- The battery must not be short circuited, as it may cause burns or become a fire hazard
- The battery may explode, catch fire or cause burns if it is incorrectly handled
- Waste batteries must be disposed of properly, in order to preserve the environment.

ENVIRONMENTAL CONDITIONS FOR USE, STORAGE AND TRANSPORT

- Do not use VEINOPLUS in the presence of working electro-surgery devices
- Do not use VEINOPLUS near (< 0.5 m) working microwave ovens
- The VEINOPLUS did not produce excessive electromagnetic interference (EMI) in tests according to EN 60601-1-2: 2001 standard. However, EMI may still occur in certain sensitive environments
- VEINOPLUS must be used in an ambient temperature range of +5 to +40°C and RH below 70%

- VEINOPLUS must be stored and transported in temperatures from -5 to +40°C and RH below 70%.

SERVICE

The VEINOPLUS device requires no routine maintenance nor calibration. When needed, for occasional cleaning / disinfecting of the plastic covers, use alcohol wipes.

All after sale services are handled by Ad Rem Technology or authorized representatives to whom, on request, the supplier will make available component lists, diagrams, descriptions and materials necessary for repairs. To obtain service, please provide the year of manufacture and a serial number of your VEINOPLUS. The year of manufacture can be read as the first four digits and the serial number as the last five digits on the label placed under the battery inside of your VEINOPLUS.

WARRANTY

VEINOPLUS is guaranteed for 12 months against any manufacturing defect, except for the battery, glass breakage and normal wear and tear of electrodes, as well as any damage due to abusive handling or removing the plastic covers of VEINOPLUS by unauthorised third parties.

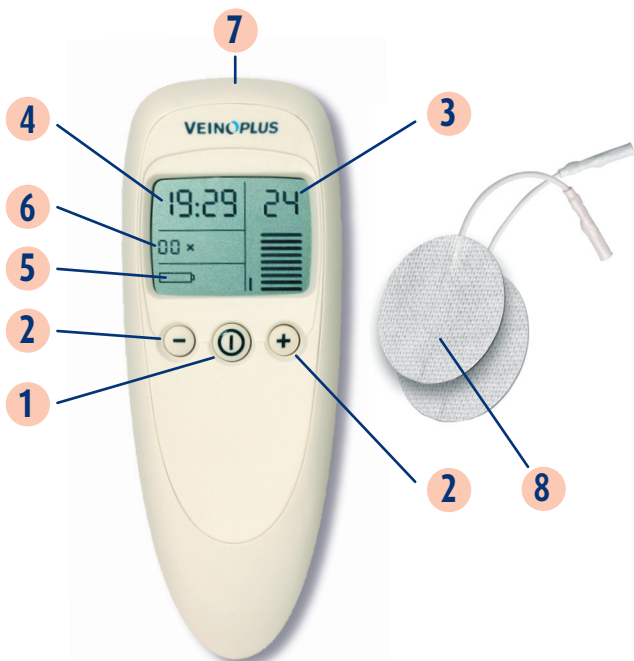
To obtain the repair or replacement of your unit under guarantee, please send your unit with its sales slip to the MANUFACTURER and EU IMPORTER:

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